

# HiPP-3 - Electron Spectrometer

## Outstanding Imaging XPS at Ambient Conditions

The HiPP-3 analyser features a newly developed technology for outstanding imaging XPS at ambient conditions with a spatial resolution of better than 10  $\mu\text{m}$ . In addition, the unique Swift Acceleration Mode enables unparalleled transmission, with count rates improved by up to a factor of 10.

Where its predecessor, the Scienta Omicron HiPP-2, was developed for HAXPES, enabling kinetic energies up to 10 keV, the Scienta Omicron HiPP-3 analyser has been optimised for XPS energies (up to 1500 eV). In this range, the transmission has been dramatically increased with the development of a new lens mode, described below.

Conceptually, Scienta Omicron HiPP-3 is similar to earlier APPES models in that it consists of a traditional Scienta Omicron R4000 analyser amended with a pre-lens that contains both differential pumping and electrostatic lenses. However, the pre-lens has been completely redesigned to offer the best transmission and spatial mode measurements at XPS energies. As before, there is an exchangeable front cone separating the pre-lens from the sample, enabling the optimum balance between pressure and transmission.

### Excellent transmission at ambient conditions

The HiPP-3 analyser features a new and improved version of the Swift Acceleration Mode. With this mode, electrons are accelerated towards the analyser by a large voltage applied to the second aperture of the differential pumping stage, see figure 1. Due to this acceleration, the inelastic scattering of the photoelectron in the surrounding gas decrease dramatically. Compared to traditional transmission modes, count rates are increased by up to one order of magnitude.

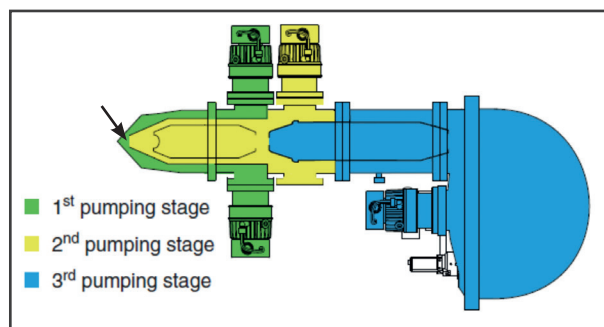


Figure 1: Differential pumping stages in Scienta Omicron HiPP-3. The arrow indicates where a high voltage is applied to accelerate electrons into analyser.

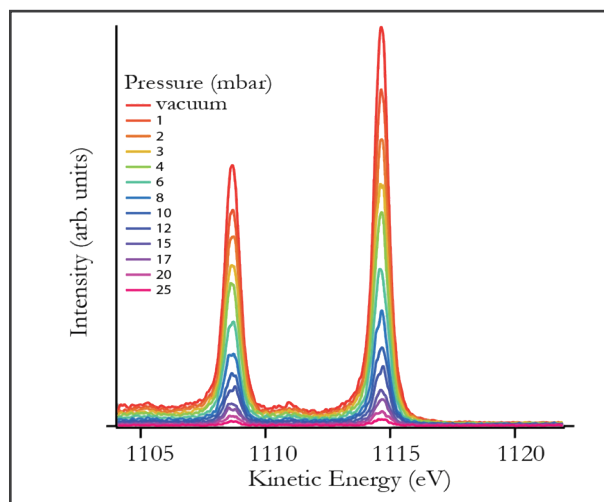


Figure 2: Ag 3d spectra recorded in different pressure of  $\text{N}_2$ . The measurement time for each spectrum was 60 seconds.

### Technical Highlights:

- Imaging XPS < 10  $\mu\text{m}$  resolution
- Swift Acceleration Mode for unprecedented transmission
- Refocusing pre-lens with efficient differential pumping
- Angular resolved range: 18°
- Interchangeable / customized front cones
- A breakthrough in Ambient Pressure PES
- Imaging XPS
- Outstanding Transmission

Transmission tests have been performed with the Scienta Omicron MX650 HP Al K $\alpha$  X-ray source. The resulting pressure dependence can be seen in figure 2. The figure shows a pressure series of Ag 3d from vacuum up to 25 mbar of N<sub>2</sub>. Due to the Swift Acceleration Mode, the signal decrease due to inelastic scattering is much less severe than for traditional analysers. An example of a spectrum recorded at 25 mbar is seen in figure 3.

This spectrum demonstrates that high resolution spectra with good statistics can be rapidly recorded even with a lab source.

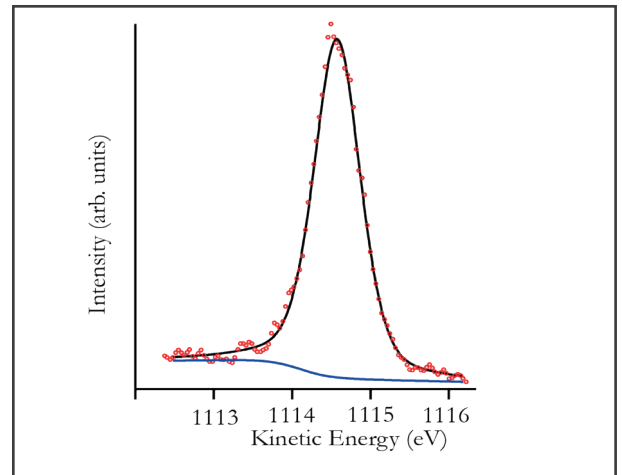


Figure 3: Ag 3d spectra recorded in 25 mbar pressure of N<sub>2</sub>. The measurement time was 17 minutes. Analyser settings corresponds to a peak width of 0.7 eV.

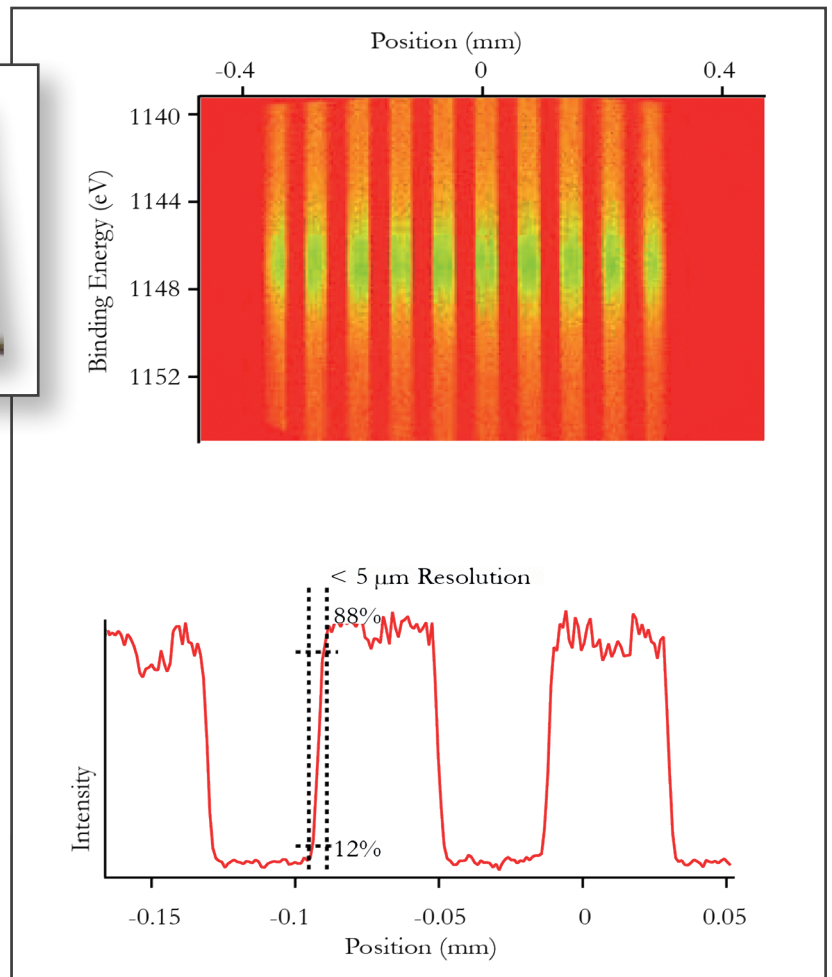
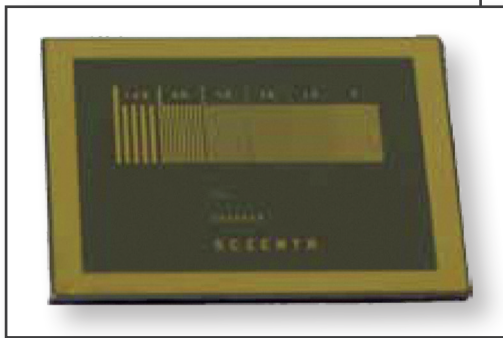



Figure 4: Demonstration of spatial resolution.  
 Top (Insert): A special test sample consisting of narrow lines of Au on Si, with the sample region magnified to the right.  
 Center: Detector image showing parallel detection of kinetic energy and emission position.  
 Bottom: Integrated spectrum showing a spatial resolution of better than 5  $\mu\text{m}$ . This resolution is achieved by a new patent pending imaging technology.

### Spatial Mode measurements

The HiPP-3 analyser features a revolutionary imaging mode. This mode enables parallel imaging in one direction with a guarantee resolution of better than 10  $\mu\text{m}$  in a field of view of up to 0.8 mm. This mode is available in both UHV and at ambient conditions.

Figure 4 show test data verifying the performance of the imaging mode. Test measurements have been performed at a dedicated development system at Scienta Omicron using the Scienta Omicron MX650 HP, which is an Al  $K\alpha$  monochromatic X-ray source adapted for high pressure, and a special test sample with narrow gold lines on silicon. In this study, gold lines with the width of 40  $\mu\text{m}$  have been used. Measurements have been performed for the Au 4d core level.

The top of figure 4 shows a picture of the test sample. The center figure show the parallel detector image, with kinetic energy along one axis and position along the other. The bottom figure show a cut out of the detector image. From the slope of the curve, a spatial resolution of better than 5  $\mu\text{m}$  has been demonstrated.

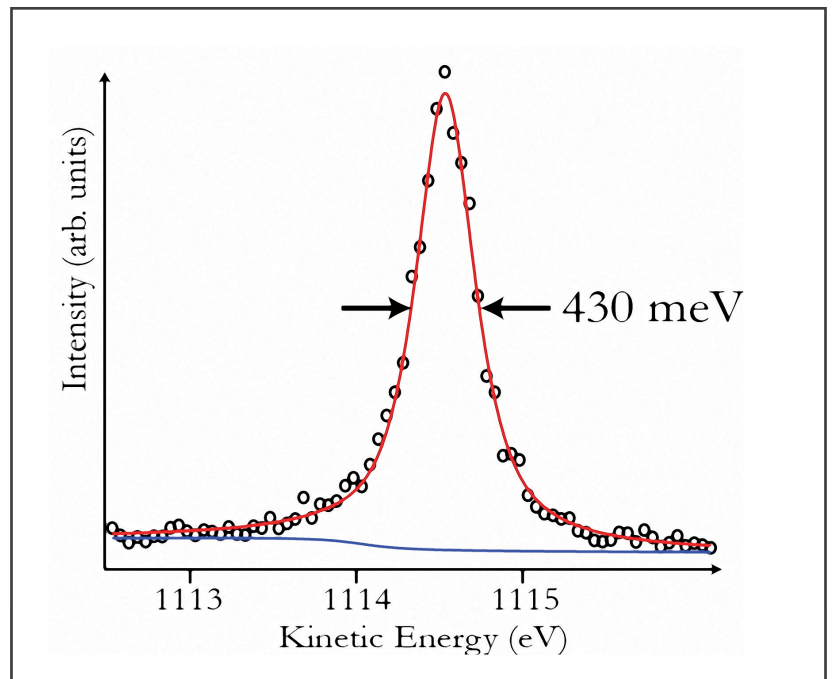

#### High Resolution XPS

With the excellent resolution and transmission of the Scienta Omicron HiPP-3 analyser in combination with the powerful monochromatic x-ray source Scienta Omicron MX650 HP, high resolution XPS can be recorded faster than ever, at both UHV and ambient conditions. The figure to the right shows a Ag 3d5/2 spectrum with a FWHM of less than 430 meV. Black circles are data points, blue is a background line and red line is the fitted spectrum.

# Technical Data

## Technology Overview

Property	Target Specification	High Voltage Electronics Property	Specification
Energy resolving power	> 1750 (0.2 mm slit)	Temperature stability	< 2 ppm/°C (R-version) (typical 0.5 ppm/°C)
Max., theoretical energy resolving power	4000 (0.1 mm slit)	Noise (AV at analyzer)	< 1 ppm + < 500 $\mu$ V (typical 0.5 ppm + < 200 $\mu$ V)
Pressure	< $5 \times 10^{-10}$ mbar Up to 30 mbar, verified for N <sub>2</sub>	Drift	< 20 ppm/year (typical 10 ppm/year)
Baking temperature	120 °C	Electric isolation	6 kV
Analyzer radius	200 mm	Min. step size HV100	1.6 mV
Mounting flange	NW 200 CF, fixed	Min. step size DAC	200 $\mu$ V
Working distance	0.3 - 9 mm, depending on aperture	DAC	Bits 16
Slits	9	Modular	Yes
Detector type	MCP/CCD camera	Communication	USB
Detector interface	Ø 40 mm MCP		
Energy channels	> 1000 simultaneous		
Angular channels	> 750 simultaneous		
Lens modes	Transmission, Swift Acceleration, Angular, Spatial		
Acquisition modes	Swept, fixed		
Detector modes	Pulsed, ADC		
Intensity deflectors	Yes (x, y)		
ISS	Option		
Analyser pumping recommendations	First stage: 2x300 l/s Second stage: 1x300 l/s Third stage: 1x300 l/s		

Property:	HiPP-3 XPS Target Specification
Energy resolution:	15 meV at 500 eV
Spatial resolution:	10 $\mu$ m meV up to 1200 eV 20 $\mu$ m meV up to 1500 eV
Angular resolved range:	$\pm 9^\circ$
Kinetic energy range	
Transmission mode:	5 - 1500 eV
Angular mode:	5 - 1500 eV
Spatial mode:	5 - 1500 eV
Pass Energy:	5 - 200 eV

Property:	HiPP-3 UPS upgrade* Target Specification
Energy resolution:	1.8 meV at 10 eV
Spatial resolution:	10 $\mu$ m meV
Angular resolved range:	$\pm 9^\circ$
Kinetic energy range	
Transmission mode:	2 - 100 eV
Angular mode:	2 - 100 eV
Spatial mode:	2 - 100 eV
Pass Energy:	2 - 10 eV

Property:	HiPP-3 6 keV upgrade* Target Specification
Energy resolution:	40 meV at 6 keV
Angular resolved range:	$\pm 9^\circ$
Kinetic energy range	
Transmission mode:	1500 - 6 keV
Angular mode:	1500 - 1950 eV
Spatial mode:	
Pass Energy:	

**How to contact us for further info:**

[www.ScientaOmicron.com](http://www.ScientaOmicron.com)  
[info@ScientaOmicron.com](mailto:info@ScientaOmicron.com)



\*Specifications in addition to HiPP-3 XPS